

# h2ofloss

## Water Flosser

### QUICK START GUIDE

H2ofloss Collapsible Reservoir Oral Irrigator




FB Group:  
h2ofloss Help & Info Hub

**Join Us Now!**

Scan to Join Our Facebook Group and Participate in Our Exciting Giveaway with a 100% Winning Rate!

Extend Your 1-Year Warranty or Get a Chance to Win a Brand New Product for Free!

## 1 Unit Prep

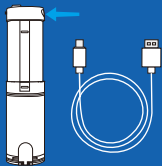


**A:** First, charging unit for 4-5 hours using a Type-C USB cable (Included)



**Note:**  
Do not press tip eject button while unit is running.

**B:** Insert tip until it clicks into place. To remove tip, pls push the tip eject button.

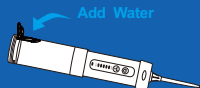

**C:** Lift the flip top on reservoir and fill with warm water or mouthwash. Refill as needed.

## 2 Water Floss



**DIY Mode Button**

**A:** Press mode button to adjust pressure. Start on the lowest pressure and increase to your liking.



**B:** Turn power ON. Slightly close lips to prevent splashing. Aim water at the gumline at a 90 degree angle



Whitening tablets not included

Button to open the nozzle

**Precision Whitening**

**C:** Upgraded multi-purpose nozzle for the H2ofloss cordless water flosser, Fill whitening tablets for enhanced oral care.

### Warning Tips:

**TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON  
PLEASE USE 5V 1A 5W ADAPTER**